
Panasonic

Portable Oral Irrigator
Irrigador oral portátil

Model No. EW1270
Modelo N.º EW1270

Operating Instructions Instrucciones de funcionamiento

Before operating this device, please read these instructions completely.
Antes de usar este dispositivo, lea completamente estas instrucciones.

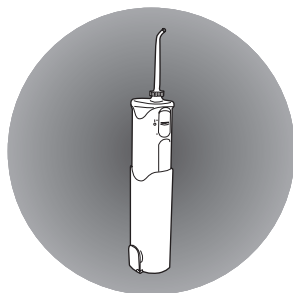Matsushita Electric Works, Ltd.
1048, Kadoma, Osaka 571-8686, Japan

No.1 EN, ES

Printed in Japan

INTRODUCTION

Thank you for purchasing the Panasonic Oral Irrigator EW1270. With its water jet of 1,600 pulses per minute, the Panasonic Oral Irrigator gently yet thoroughly cleanses your teeth and gums, removing trapped food particles and reducing bacteria that can lead to gingivitis and other periodontal diseases. When used in combination with regular brushing and flossing, the Panasonic Oral Irrigator can help you have healthier teeth and gums.

Before using this device, please read these instructions completely. Keep these instructions for reference.

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IMPORTANT SAFETY PRECAUTIONS

1. Do not use while bathing.
2. Do not place the device where it can slide, fall or be pulled into water.
3. Close supervision is necessary when this device is used by, on or near children.
4. KEEP JET NOZZLE ATTACHMENT AWAY FROM SMALL CHILDREN. The nozzle is small enough to be swallowed and can become lodged in the throat.
5. Keep this device out of the reach of children when not in use.
6. This device is designed for use only in the oral cavity. Do not direct water in delicate areas, such as under the tongue, in the eye, ear or nose, or in other delicate areas of the body.
7. Use only attachments recommended by the manufacturer.
8. DO NOT operate this product if it is not working properly or it has been dropped or damaged. If damaged or inoperative, return the product to a service center for inspection and repair.
9. NEVER drop or insert any object into any opening on the device.
10. Do not disassemble the device. Damage to the product and/or personal injury may occur.
11. People with tooth or gum disease or people receiving oral treatment should consult their dentist before use.
12. People who experience tooth or gum pain may not be able to use this device and should consult their dentist.
13. NEVER use while drowsy.
14. Remove the batteries from the device if it will not be used for more than a month. Otherwise, electrical leakage may occur if the insulator is worn.
15. Fill the tank with water only.
16. Certain types of mouthwash may damage the equipment by causing the housing or water tank to crack. Therefore, do not use mouthwash in the water tank.
17. Improper handling of batteries may result in battery explosion or in corrosion from battery leakage. Please observe the following to ensure proper use of batteries.
 - a. Be sure to turn off the power of the device after use.
 - b. Do not mix different types or sizes of batteries.
 - c. Replace all batteries at the same time. Do not mix old and new batteries.
 - d. Be sure to insert batteries with correct polarity, as instructed.
 - e. Remove batteries when they are worn out, and dispose of them properly according to all applicable environmental regulations.
 - f. Do not disassemble or throw batteries into a fire.
 - g. Do not short-circuit batteries.
 - h. Do not charge the batteries included with this unit.

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BEFORE USE

Loading Batteries

• Be sure to wipe the unit before installing the batteries so that no water enters the battery compartment.

1. Twist the battery cover to the OPEN position and open it (Fig. 2).
2. Insert two 1.5-V alkaline (AA/LR6) batteries, ensuring correct "+" and "-" polarity.

• The device can be used approximately 30 times (with one full tank of water per use) with a fresh set of Panasonic alkaline batteries.

• Do not mix old and new batteries.

3. Twist the battery cover to the CLOSE position and close it.

Note:

The batteries' power output is impaired at low temperature.

CAUTION:

Do not throw batteries into a fire. Doing so may result in an accident or explosion.

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USING YOUR ORAL IRRIGATOR

(Figs. 3 & 4)



• Rinsing your mouth with mouthwash is not a substitute for proper dental care. Brush your teeth or dental prosthesis for at least five minutes immediately after every meal, then use the device.

1. Open the water tank lid.
 - The tank cannot be pulled out easily when the lid is closed.
2. Pull the water tank section downward to the water tank usable line.
3. Remove the nozzle from its compartment and insert firmly into the nozzle insertion hole until the nozzle clicks into place.
4. Fill the tank with cold or warm water, not warmer than 40°C (104°F).
 - Do not put any toothpaste, powder or other chemicals into the water tank.
 - Certain types of mouthwash may damage the unit by causing the housing or tank to crack. Therefore, use only water in the water tank.
5. Close the lid.

CAUTION:

NEVER fill the water tank with hot water! Sensitive gums can be easily irritated or damaged by hot water.

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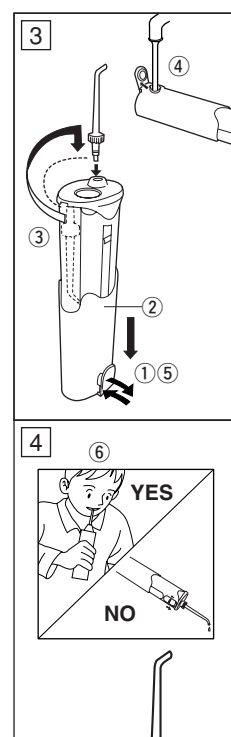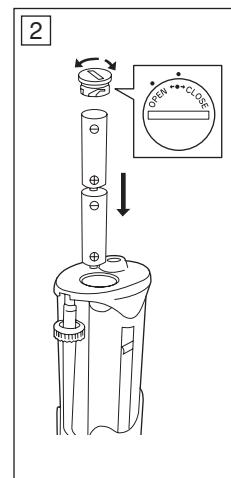PARTS IDENTIFICATION (Fig. 1)

Main unit

- ① Battery cover
- ② Nozzle insertion hole
- ③ Switch
- ④ Nozzle compartment
- ⑤ Water tank

• There may be moisture inside the water tank and suction hose. This is moisture left from the distilled water used for product testing at the factory and is hygienically safe.

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
Accesorio

- ⑦ Boquilla (1)

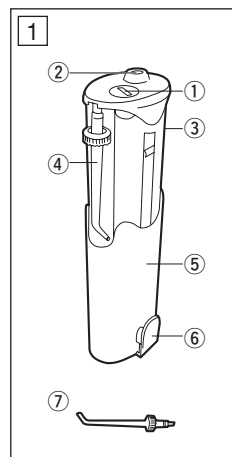6. Lean over the sink and place the nozzle in your mouth. Close your lips slightly to prevent splattering.

- Hold the device upright when in use. Excessive tilting can prevent the water from coming out of the nozzle.
- The connecting section of the main unit and the water tank is not waterproof. When the tank is filled with water, be sure to place the device upright. Also, be sure to set the device down gently so that the bottom of the tank is not subjected to strong impact, which can cause the lid to open and the water to spill out.

7. Slide the switch up to turn on the device.

PRECAUTIONS:

- Do not turn the switch on before putting the nozzle in your mouth.
- Do not press the tip of the nozzle too firmly against your teeth or gums. Doing so may injure your gums.
- If you are suffering from tooth or gum disease or receiving oral treatment, consult a dentist before use.
- The device is designed exclusively for use in the oral cavity. Never use it for any other purpose.
- During the first several uses, light bleeding of the gums may occur. If your gums are healthy, any bleeding should stop in one to two weeks. If bleeding persists for two or more weeks, consult your dentist before continuing use.

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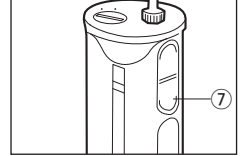
TIPS ON USE

- Bend over the sink and hold the unit in an upright position with the nozzle tip pointing toward your teeth.
- Keep your mouth opened slightly so that the water can flow out.
- Starting in the back (molar) area and working toward the front teeth, systematically rinse all teeth and interdental spaces by:
 - directing the water jet flow at a right angle to the teeth and the area between the teeth and gums;
 - directing the water jet slowly along the teeth;
 - directing the water jet along the gumline with the nozzle tip at a right angle to the gum tissue;
 - cleaning between the teeth and around any braces, crowns or bridgework thoroughly.

Notes:

One full tank of water is not enough to rinse all of your teeth and gums adequately. We therefore recommend refilling the tank several times to ensure thorough cleaning of teeth and gums.

- In addition to using the oral irrigator and brushing and flossing your teeth, we recommend regular checkups by a dentist as part of a complete dental hygiene program.

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AFTER USE (Fig. 5)

- Open the water tank lid and pour out the water.


Caution: NEVER immerse the device in water.

- Turn on the device for a few seconds until any water remaining inside the device and nozzle drains out completely.

Note:

A small amount of water may remain inside the tank even after you have drained it by turning on the device for a few seconds. Be aware of such possible water drainage when you put away the device or pack it when traveling.

- Remove the nozzle and place it into the nozzle compartment.
- Push the tank up.
 - Be careful not to accidentally insert or pinch your finger.
 - When pushing the tank up, leave the lid open and make sure that the tip of the nozzle is inside the tank.
- Close the water tank lid.
 - After rinsing the unit, wipe and dry it well before storing the unit to prevent contamination and bacterial growth.

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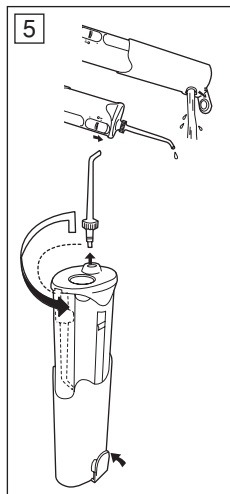TROUBLESHOOTING – Possible causes and Remedies

Problem	Possible cause	Remedy
The device stops immediately after being switched on or does not operate at all.	Batteries are not inserted properly.	Insert the batteries properly. See the Loading Batteries section.
	Battery power is low.	Replace the batteries with new ones. See the Loading Batteries section.
Water pressure is too low or operating time is very short.	Battery power is low.	Replace the batteries with new ones. See the Loading Batteries section.
The pump is running but no water comes out of the nozzle.	The tank is empty.	Fill up tank with water.
	The tank is not pulled down enough.	Pull down the tank fully.
	Wrong position.	Use the device in an upright position with the nozzle pointing upwards.
	The nozzle is clogged or deformed.	Replace the nozzle.

If the device still does not function properly, stop using and contact an authorized service center or the store from which the device was purchased.

CAUTION:

- Never attempt to repair the device yourself.
- If parts must be exchanged or replaced, or if the device malfunctions, contact an authorized service center.

SPECIFICATIONS

Power supply: 3 V DC (Two 1.5-V alkaline batteries (AA/LR6))

Pump frequency: 1600 pulses/minute

Water pressure: approx. 440 kPa (approx. 4.5 kgf/cm²)

Battery life: approx. 20 minutes (with Panasonic alkaline AA/LR6 batteries)

Reservoir volume: 165 ml (approx. 40 seconds of use with new batteries)

Dimensions: 225(H) x 45(W) x 63(D) mm (9"(H) x 1-3/4"(W) x 2-1/2"(D))
150(H) x 45(W) x 63(D) mm (6"(H) x 1-3/4"(W) x 2-1/2"(D)) when the tank is closed.

Weight: 200 g (7 oz) excluding batteries

Specifications are subject to change without notice.

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MAINTENANCE

① Nozzle

Remove the nozzle from the device. Rinse with water and wipe with a dry cloth.

② Main unit

Clean the main unit by wiping it with a cloth lightly dampened with liquid soap.

- Use only water or a neutral cleaning agent to clean the main unit. Never use caustic or abrasive agents (e.g., vinegar or lime-scale remover) because they may damage the device.

③ Tank


Pull the tank down firmly to remove it. Clean the inside with a soft brush or similar object if necessary. Wipe off any residue inside the tank. If the device will not be used for more than a week, be sure to dry the inside of the tank.

④ Water tank lid

Wipe off to prevent water leakage.

⑤ Hose

Do not pull or bend the hose forcibly. Doing so may damage the hose and its spring.

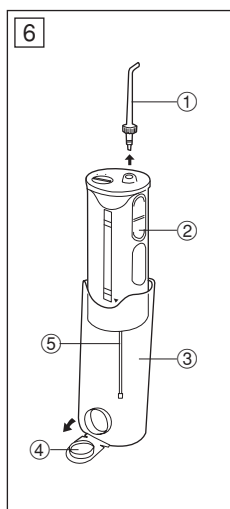CAUTION:

- Do not wash with hot water (exceeding 50°C or 122°F).
- Do not immerse the main unit in water.

Note:

If the device will not be used for more than a month, be sure to remove the batteries.

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TRAVELING WITH THE PANASONIC ORAL IRRIGATOR

Thanks to its handy shape and compact size, the Panasonic Oral Irrigator is an ideal traveling companion. When using it on a trip, be sure to completely drain the tank and nozzle before putting the device away.

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